



A Dräger and Siemens Company

Field Service Procedure

Part Number: SP00282

Rev: -

Date: 23 July 2003

© 2003 Draeger Medical, Inc.

Patient Airway Gauge Restrictor Installation on Fabius GS®

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)
[RETURN TO CD-ROM TABLE OF CONTENTS](#)

Patient Airway Gauge Restrictor Installation on Fabius GS

INSTALLATION PROCEDURE

1. Verify unit is in 'Standby' mode and all flow control valves are closed.
2. At the tee fitting, disconnect the patient airway (gauge) hose. Refer to Figure 1.
3. Insert restrictor into patient airway gauge hose approximately one (1) inch from end of tube and align the restrictor orifice with the airway.
4. Reconnect patient airway gauge hose to tee fitting.
5. Perform Breathing System Leak test in accordance with instructions provided in paragraph 7.9.3 of the Fabius GS Technical Service Manual.

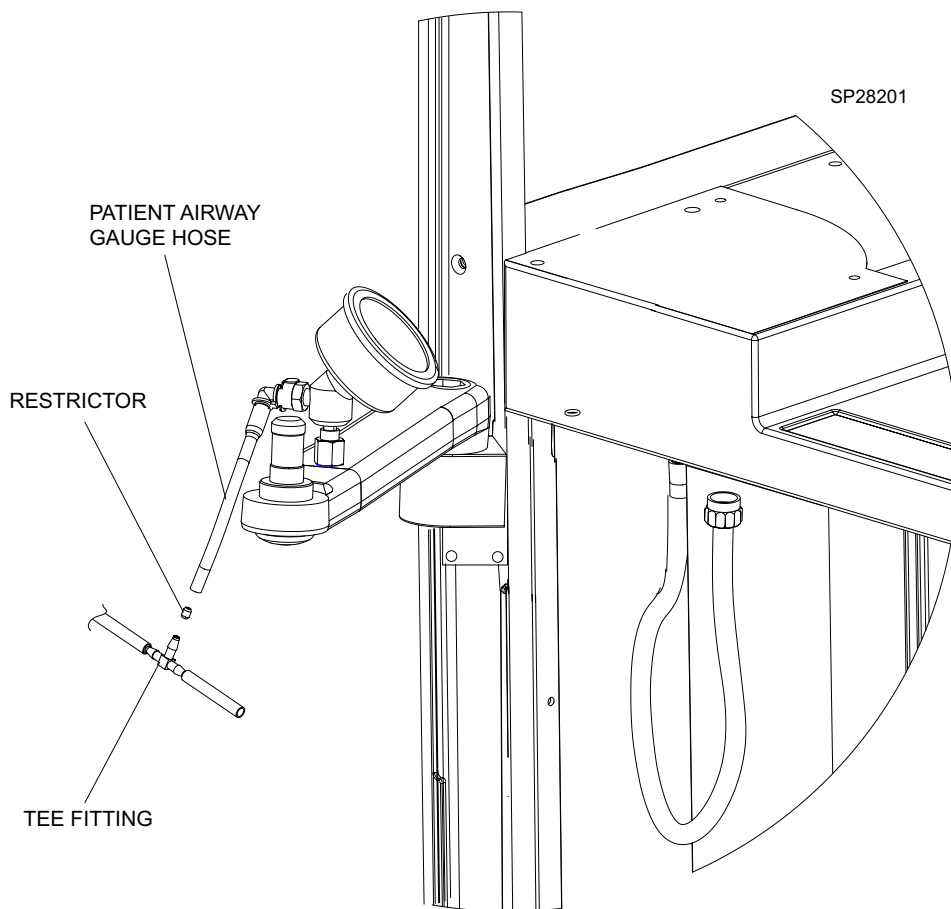


Figure 1. Restrictor Installation

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)
[RETURN TO CD-ROM TABLE OF CONTENTS](#)

[RETURN TO SERVICE PROCEDURE TABLE OF CONTENTS](#)
[RETURN TO CD-ROM TABLE OF CONTENTS](#)



DrägerService is a division of
Draeger Medical, Inc.
3122 Commerce Drive
Telford, PA 18969
Tel: (215) 721-5402
(800) 543-5047
Fax: (215) 721-5784
Web: www.draegermedical.com
Printed in the U.S.A.